

FIELD CHANGE ORDER

Service

Issued by : PMC Hamburg
Released : D. Morgenstern

Ref. No. : 00 186 003
Date : March 1994
Product group : 740



H.T. cables with connectors O3 / R3 / R5 / R10

APPLIES TO:

High-tension cables with plugs O3 / R3 / R5 / R10

Type numbers:	9806 400 2..02	R3 / R3
	9806 401 2..02	O3 / R5
	9806 401 7..02	O3 / R3
	9806 402 6..02	O3 / O3
	9806 402 8..02	O3 / R5
	9806 405 0..02	R10 / R10

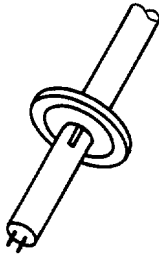
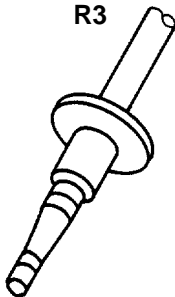
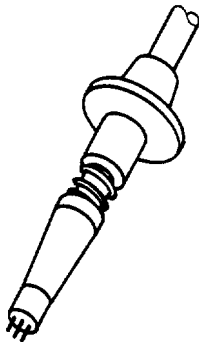
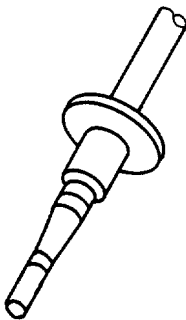
TITLE:

H.T. connectors O3 / R3 / R5 / R10

LIST OF PAGES AND DRAWINGS:

1	(94.0)
2	(94.0)

To avoid any problems during installation or an exchange of the tube refer to the table below. This table explains how to use silicone paste and silicone washers at the tube and the generator end for high-tension cables with plugs O3/R3/R5/R10.

HT plug type	HT generator	X-ray tube
<p>O3</p>  <p>Note Do not use the O-ring for this plug.</p>	<p>Silicone washer NO</p> <p>Silicone paste NO</p> <p>Fill only the HT receptacle about 1 cm with HT generator oil.</p> <p>Note With older (modular) generators provide the front third of the plug with an even, thin coat of silicone paste.</p>	<p>Silicone washer YES</p> <p>Silicone paste YES</p> <p>Bend apart the contact pins before connecting the plug. Clean the plug with alcohol only. Do not grease the HT plug with silicone paste. Grease the silicone washer with silicone paste. Do not grease the contacts.</p> <p>Note Exchange the silicone washer when exchanging the tube.</p>
<p>R3</p> 	<p>Silicone paste YES</p> <p>Clean the plug with silicone paste only. Do not grease the contacts. Provide the clean rubber cone with an even, thin coat of silicone paste.</p>	<p>Silicone paste YES</p> <p>Clean the plug with silicone paste only. Do not grease the contacts. Provide the clean rubber cone with an even, thin coat of silicone paste.</p>
<p>R5</p> 	<p>Not present at the generator end.</p>	<p>Silicone paste YES</p> <p>Clean the plug with silicone paste only. Do not grease the contacts. Provide the clean rubber cone with an even, thin coat of silicone paste.</p>
<p>R10</p> 	<p>Silicone paste YES</p> <p>Clean the plug with silicone paste only. Do not grease the contacts. Provide the clean rubber cone with an even, thin coat of silicone paste.</p>	<p>Silicone paste YES</p> <p>Clean the plug with silicone paste only. Do not grease the contacts. Provide the clean rubber cone with an even, thin coat of silicone paste.</p>